



NEW RIGID CORSETS
PRS666 / PRS660 / PRS661


Prim spine

FOCUSBRACE

*New bivalve-shaped LSO and TLSO orthoses.
Easy to fit without the need for moulds.*

Designed to limit mobility while maintaining
comfort to improve treatment adherence.

Immobilization and comfort, in harmony.

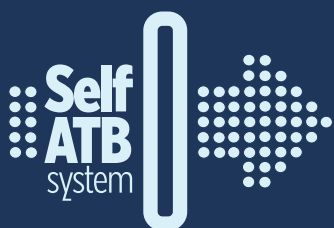

Stabilization level

Primspine FOCUSBRACE

NEW RIGID CORSETS

The intra-abdominal pressure exerted by the structure of the **PRIMSPINE Focusbrace** corrects hyperlordosis by achieving a high degree of lumbar stabilisation, immobilisation and unloading.

Bivalve-shaped braces

The posterior structure overlaps laterally on the anterior structure in a safe way, through a micro-hook.
Easy to place, even in post-surgical and bedridden patients.



Front plate

Specially designed to fit any shape abdomen.


Front structure

Anteroposterior adjustment straps

Front-to-back adjustment elastic straps that join its exclusive **Self-ATB System** buckle system.



Rear structure

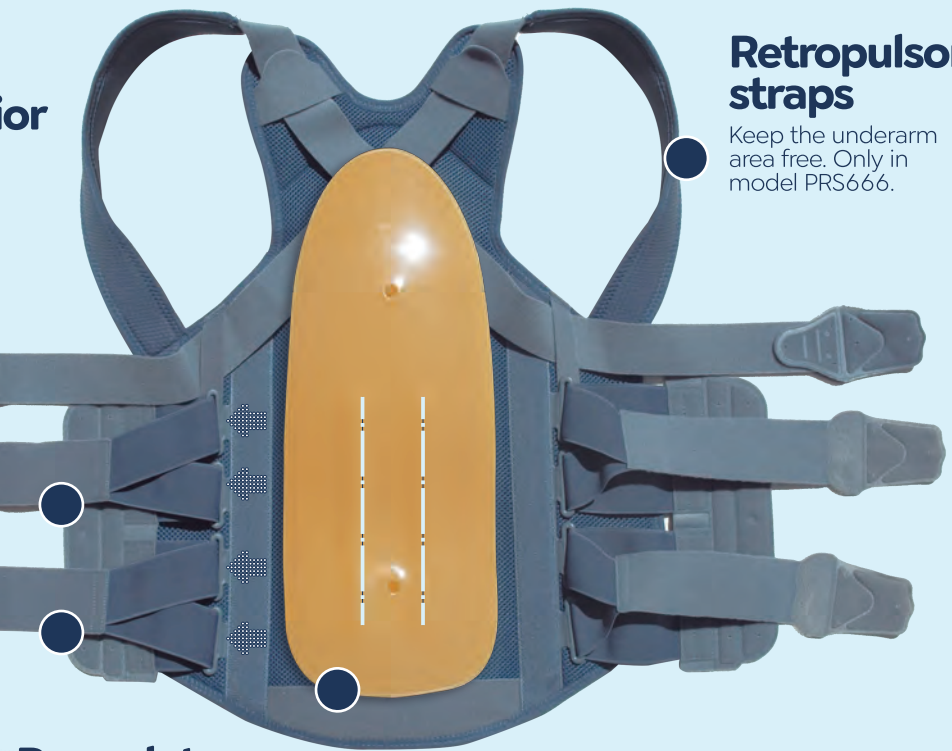


Rear plate

With paravertebral reinforcement ribs positioned between 2 metal stays for improved lumbar stabilisation and immobilisation. Mouldable with a hot air gun.

Retropulsor straps


Keep the underarm area free. Only in model PRS666.

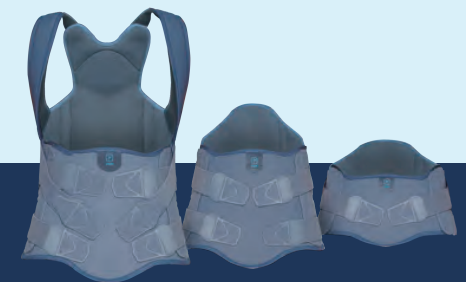

Immobilization and comfort, in harmony.

Prim spine

FOCUSBRACE

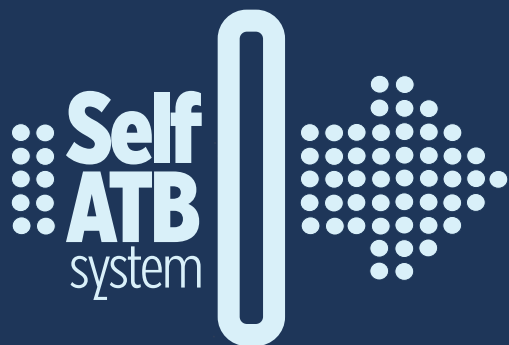

Buckles system that allows the uniform distribution of tension, favouring fitting and comfort.

3D padded fabric

Strong, breathable 3D padded fabric, with cushioning and comfortable cotton fabric in contact with the skin.



Breathable
Latex-free





Primospine FOCUSBRACE


Rigid Lumbar Sacral Orthosis (LSO) High




Size/Size range	A cm.	B cm.
S 80-90	27,5	39
M 88-100	27,5	39
L 98-110	28,5	39
XL 108-120	29,5	39

Indications:

Spinal disc herniation, osteoporosis, spondyloarthritis, spondylolysis (with or without stability), trauma, compression fractures, degenerative disc disease, postoperative fusion and chronic muscle disease.

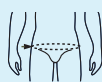

Measurement: circumference at height of pelvis

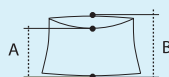

A: Front height
B: Back height





Primospine

FOCUSBRACE


Rigid Lumbar Sacral Orthosis (LSO)
Low




PRS661 ●

Size/Size range	A cm.	B cm.
S 80-90	19	28
M 88-100	19	28
L 98-110	19	28
XL 108-120	19	28

Indications:
Spinal disc herniation, osteoporosis, spondyloarthritis, spondylolysis (with or without stability), trauma, compression fractures, degenerative disc disease, postoperative fusion and chronic muscle disease.

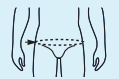

Measurement: circumference at height of pelvis

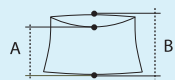

A: Front height
B: Back height



Primspine

FOCUSBRACE

Rigid Thoracic Lumbar Sacral Orthosis (TLSO)



PRS666


Size/Size range	A cm.	B cm.
S 80-90	27,5	51
M 88-100	27,5	51
L 98-110	28,5	51
XL 108-120	29,5	51

Indications:

Spinal disc herniation, osteoporosis, spondyloarthrosis, spondylolysis (with or without stability), trauma, compression fractures, degenerative disc disease, postoperative fusion, chronic muscle disease, and prevention of canal stenosis in bedridden patients.



Measurement: circumference at height of pelvis

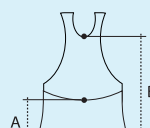

A: Front height
B: Back height

